Azura

► **Column Thermostat CT 2.1**
User Manual

V6810

HPLC

Table of Contents

Note: For your own safety, **read** the manual and **always** observe the warnings and safety information on the device and in the manual!

Intended Use	5
Front View	5
Rear View	6
Features	6
Scope of Delivery	7
Safety for Users	7
Signal Words	8
Decontamination	9
Symbols and Signs	9
Installation	10
Transport	10
Operating Environment	11
Space Requirements	11
Setup	11
Unpacking	11
Power Supply and Connection	12
Connectors for Control	12
Connecting the Device to the Computer	13
Configuring the LAN Settings	13
Connecting the Cables	14
Configuring the Router	14
Integrating the LAN into a Company Network	15
Controlling Several Systems Separately in a LAN	15
Changing the Swing Direction of the Door	16
Installing the Cartridge for Precolumn Tempering	16
Installing the Cartridge of the Precolumn Tempering	16
Operation	17
Initial Start-Up	17
Switch-On	17
Controlling with Chromatography Software	18

Controlling with Operating Panel	18
Self-Testing	18
Remote Operation Using Chromatography Software	19
Leak Sensor Settings	19
Setting the T-Lock Temperature	19
Entering a Fixed Temperature	20
Changing the Temperature	21
Entering Timed Temperature Programs	22
Programming a Temperature Gradient	23
Starting a Program or a Sequence	25
Functionality Tests	26
Troubleshooting	26
LAN	26
Possible Problems and Rectifications	27
Error messages	28
Maintenance and Care	28
Contact the Technical Support	29
Maintenance Contract	29
Cleaning and Care	29
Storage	29
Disconnecting from Power Supply	29
Removing a Leak	30
Correcting the Temperature	30
Temperature - Auto-Calibration	30
Correcting the Temperature Manually	31
Adjusting the Leak Sensor	31
Technical Data	32
Repeat Orders	34
Modules	34
Accessories and Spare Parts	34
Disposal	34
Legal Information	35
Warranty Conditions	35
Transport Damage	35
HPLC Glossary	36
Index	38
Declaration of Conformity	40

Intended Use

Note: Only use the device for applications that fall within the range of the intended use. Otherwise, the protective and safety equipment of the device could fail.

Description The Column Thermostat CT 2.1 is suitable for use with up to five columns with a maximum length of 350 mm and a maximum diameter of 20 mm.

Scope of application The device can be used in the following areas:

- Biochemistry analysis
- Chemical analysis
- Food analysis
- Pharmaceutical analysis
- Environmental analysis

Place of installation Always use the device in rooms that are well-ventilated. Do not place the column thermostat near heat sources or in direct sunlight.

Front View

Legend

- ① Pre-column tempering
- ② Door
- ③ Housing
- ④ Column holder
- ⑤ Fan
- ⑥ Temperature sensor
- ⑦ Column
- ⑧ Leak sensor

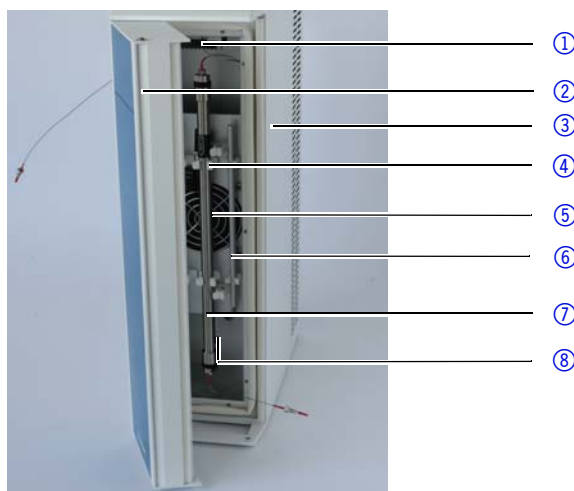
Fig. 1 Front view

Rear View

Note: On the rear side of the thermostat are situated:

- LAN connector
- Optional operating panel
- Power plug
- Serial number of the device
- Power switch (on/off)
- Potentiometer for leak sensor settings

Service interface

Please note that the service interface is solely used for repair and maintenance tasks performed by service technicians.

Legend

- ① Connector operating panel
- ② Potentiometer
- ③ LAN port
- ④ Service interface
- ⑤ Fan
- ⑥ Power switch
- ⑦ Label with serial number
- ⑧ Power plug

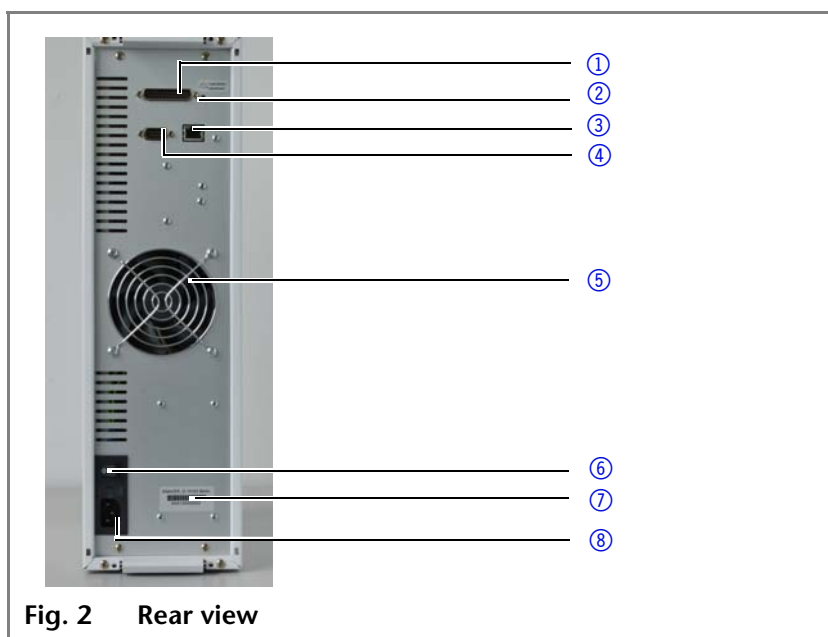
Fig. 2 Rear view

Features

The column thermostat is used to regulate the temperature of one or more HPLC columns. For this purpose, you can select a constant temperature ranging from 5 to 85 °C. As an additional option, an element for pre-column tempering is available. With this element, the temperature of the mobile phase can be adjusted to the set temperature before it enters the column.

GLP data

The different software products (e.g. ClarityChrom® or OpenLAB®) can be used to display or read out the GLP data of the thermostat. Find a detailed description on how to display or read out GLP data in the respective user manuals.

Optional accessories

The following optional accessories are available:

- Operating panel for the thermostat
- Pre-column tempering element

Scope of Delivery

Note: Only use original spare parts and accessories made by KNAUER or a company authorized by KNAUER.

Delivery	Column thermostat CT 2.1	<input type="checkbox"/>
	Accessory kit CT 2.1	<input type="checkbox"/>
	Installation Qualification	<input type="checkbox"/>

Safety for Users

Professional Group The user manual addresses persons who are qualified as chemical laboratory technicians or have completed comparable vocational training.

The following knowledge is required:

- Fundamental knowledge of liquid chromatography
- Knowledge regarding substances that are suitable only to a limited extent for use in liquid chromatography
- Knowledge regarding the health risks of chemicals
- Participation during an installation of a device or a training by the company KNAUER or an authorized company.

If you do not belong to this or a comparable professional group, you may not perform the work described in this user manual under any circumstances. In this case, please contact your superior.

Safety Equipment When working with the device, take measures according to lab regulations and wear protective clothing:

- Safety glasses with side protection
- Protective gloves
- Lab coat

What must be taken into account?

- All safety instructions in the user manual
- The environmental, installation, and connection specifications in the user manual
- National and international regulations pertaining to laboratory work
- Original spare parts, tools, and solvents made or recommended by KNAUER
- Good Laboratory Practice (GLP)
- Accident prevention regulations published by the accident insurance companies for laboratory work
- Filtration of substances under analysis
- Use of inline filters

- Once they have been used, never re-use capillaries in other areas of the HPLC system.
- Only use a given PEEK fitting for one specific port and never re-use it for other ports. Always install new PEEK fittings on each separate port.
- Follow KNAUER or manufacturer's instructions on caring for the columns

More safety-relevant information is listed below:

- flammability: Organic solvents are highly flammable. Since capillaries can detach from their screw fittings and allow solvent to escape, it is prohibited to have any open flames near the analytical system.
- solvent tray: Risk of electrical shock or short circuit if liquids get into the device's interior. For this reason, place all bottles in a solvent tray.
- solvent lines: Install capillaries and tubing in such a way that liquids cannot get into the interior in case of a leak.
- leaks: Regularly check if any system components are leaking.
- power cable: Defective power cables are not to be used to connect the device and the power supply system.
- self-ignition point: Only use eluents that have a self-ignition point higher than 150 °C under normal ambient conditions.
- power strip: If several devices are connected to one power strip, always consider the maximum power consumption of each device.
- power supply: Only connect devices to voltage sources, whose voltage equals the device's voltage.
- toxicity: Organic eluents are toxic above a certain concentration. Ensure that work areas are always well-ventilated! Wear protective gloves and safety glasses when working on the device!

Where is use of the device prohibited?

Never use the system in potentially explosive atmospheres without appropriate protective equipment. For further information, contact the Technical Support of KNAUER.

Decommissioning the Device Securely

At any time, take the device completely out of operation by either switching off the power switch or by pulling the power plug.

Opening the Device

The device may be opened by the KNAUER Technical Support or any company authorized by KNAUER only.

Signal Words

Possible dangers related to the device are divided into personal and material damage in this user manual.

DANGER! Lethal or very serious injuries can occur.

WARNING! Serious or moderate injuries can occur.

CAUTION! Moderate injuries and/or device defects can occur.

Decontamination

Contamination of devices with toxic, infectious or radioactive substances poses a hazard for all persons during operation, repair, sale, and disposal of a device.

DANGER!

Health danger if getting in contact with toxic, infectious or radio-active substances.

Before disposing of the device or sending it away for repair, you are required to decontaminate the device adequately.

All contaminated devices must be properly decontaminated by a specialist company or the operating company before they can be recommissioned, repaired, sold, or disposed of. All materials or fluids used for decontamination must be collected separately and disposed of properly.

Symbols and Signs

The following symbols and signs can be found on the device, in the chromatography software or in the user manual:

Warning signs

Symbol	Meaning
	Electric shock hazard
	General warning sign, moderate injuries can occur and also damages to system, device, or components.
	Beware of explosion due to high concentration of organic solvent vapor.
	Electrostatic discharge hazard, damages to system, device, or components may occur.
	Possible damage caused by leaks at the thermostat.
	Hazard caused by open thermostat door.

CE mark

Symbol	Meaning
	Hazardous substances.
	Hot-surface hazard.
	A device or system marked with CE fulfills the product specific requirements of European directives. This is confirmed in a Declaration of Conformity.
	Testing seals in Canada and the USA at nationally recognized testing centers (NRTL). The certified device or system has successfully passed the quality and security tests.

Installation

The chapter Installation describes all preparatory steps prior to the start-up. If you encounter difficulties during installation, contact the Technical Support.

Contact data

Phone	+49 30 809727-111
Fax	+49 30 8015010
E-mail	support@knauer.net

Transport

Carefully prepare the device for transport or storage. If you want to return your device to KNAUER for repairs, enclose the Service Request Form which can be downloaded from our website.

Device data

For a secure transport, note the weight and dimensions of the thermostat (see Technical Data).

CAUTION!

Damage to the device by carrying or lifting it on protruding housing parts.

Lift the device on the side of the housing only.

Lifting

Clasp the device at its side panels and lift it out of the packaging. Do not hold onto front cover.

Operating Environment

Ambient Conditions Only if the requirements for ambient conditions are met, can the intended use be ensured.

CAUTION!

Defect of the device due to overheating possible.

- Set up the device in such a way that it is protected against exposure to direct sunlight.
- Keep at least 15 cm clear at the rear and 5–10 cm at each side for air circulation.

- Sunlight: Protect the device against direct exposure to sunlight.
- AC system: Set up the device at a location not exposed to air drafts.
- Installation site: Position the device on a level and even surface.
- Humidity < 90 %, non-condensing
- Temperature 10-40 °C (50-104 °F)
- Recommended working environment 17-28 °C (62.6-82.4 °F)
- Height above sea level maximum 2000 meters

Space Requirements

To protect the device from overheating, allow enough space on each side of the device:

- At least 5 cm if there is another device on one side.
- At least 10 cm if other devices set up on both sides.
- At least 15 cm on the rear panel for the fan.

Setup

After you have selected a suitable place for operation, remove the device packaging. The requirements and a description can be found in the following section.

Unpacking

Prerequisite Check packaging for damage caused during transportation.
Tool Utility knife

CAUTION!

Damage to the device by carrying or lifting it on protruding housing parts.

Lift the device on the side of the housing only.

Procedure

Process

1. Set-up the package in such a way that you can read the label. Using the utility knife, cut the adhesive tape.
2. Remove the foam insert. Take out the accessories kit and the manual.
3. Open the accessories kit and check the scope of delivery. In case any parts are missing, contact the Technical Support.
4. Clasp the device at its bottom and lift it out of the packaging. Do not hold onto front cover.
5. Check the device for signs of damage that occurred during transport. In case you notice any damage, contact the Technical Support.
6. Set-up the device in its location.

Next steps

Store packaging and keep the included packing list for repeat orders.

Power Supply and Connection

The device is suitable only for a voltage range of 90-230 Volt AC and 50-60 Hz. The devices are switched on and off with the power switch on the rear side.

Prerequisites

- The electrical power supply at the installation site must be connected directly to the nearest main power line.
- The power must be free from ripple, residual current, voltage peaks and electromagnetic interference.
- The connectors for the mains voltage are grounded accordingly.
- The device receives sufficient power with reserve capacity.

Power cable

Use only the enclosed power cable to connect the device to the mains. Replace defective power cables only with accessories from KNAUER. Only use power cables admitted for use in your country.

Power plug

Make sure that the power plug on the rear of the device is always accessible, so that the device can be disconnected from the power supply.

Connectors for Control

The detector can be controlled in two ways:

- Via an operating panel (optional)
- As part of a LAN, via the LAN connector of the router

All connectors for external control are located on the rear side of the thermostat.

With an operating panel, the thermostat can be controlled locally. The operating panel has to be bought separately. The operating panel is connected via a serial interface on the back of the device.

Legend

- ① Connector operating panel
- ② LAN connector
- ③ Service interface

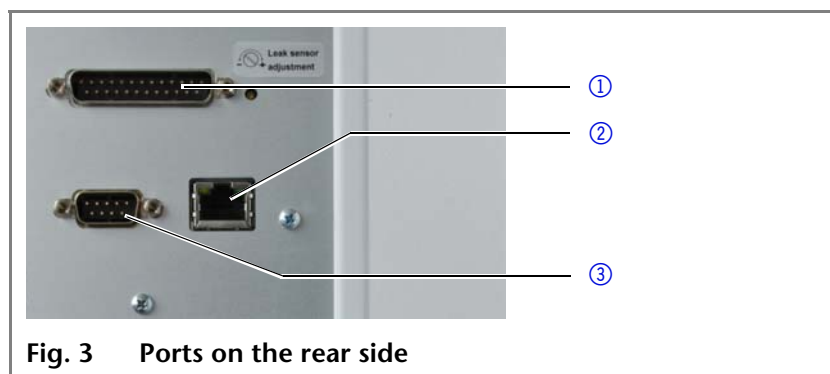
Fig. 3 Ports on the rear side

Note: While the thermostat is switched on, do not remove or connect the operating panel. Always switch off the thermostat to do so. Fixate the plug of the operating panel with the knurled-head screws.

Connecting the Device to the Computer

This section describes how to set up an HPLC system in a local area network (LAN) and how a network administrator can integrate this LAN into your company network. The description applies to the operating system Windows® and all conventional routers.

Note: To set up a LAN, we recommend to use a router. That means the following steps are required:

- Process**
1. On the computer, go to the control panel and check the LAN properties.
 2. Hook up the router to the devices and the computer.
 3. On the computer, configure the router to set up the network.
 4. Install the chromatography software from the data storage device.
 5. Switch on the device and run the chromatography software.

Configuring the LAN Settings

The LAN uses only one server (which is normally the router) from that the devices automatically receive their IP address.


- Prerequisite**
- In Windows®, power saving, hibernation, standby, and screen saver must be deactivated.
 - In case you use an USB-to-COM box, the option "Allow the computer to turn off this device to save power" in the device-manager must be deactivated for all USB hosts.

- Only for Windows 7: For the network adapter, the option "Allow the computer to turn off this device to save power" in the Device Manager must be deactivated.

- Procedure**
1. In Windows 7 choose *Start* ⇒ *Control Panel* ⇒ *Network and Sharing Center*.
 2. Double-click on *LAN Connection*.
 3. Click on the button *Properties*.
 4. Select *Internet Protocol version 4 (TCP/IPv4)*.
 5. Click on the button *Properties*.
 6. Check the settings in the tab *General*. The correct settings for the DHCP client are:
 - a) *Obtain IP address automatically*
 - b) *Obtain DNS server address automatically*
 7. Click on the button *OK*.

Connecting the Cables

A router ③ has several LAN ports ② and one WAN port ④ that can be used to integrate the LAN into a wide area network (WAN), e.g. a company network or the Internet. In contrast, the LAN ports serve to set up a network from devices ① and a computer ⑤. To avoid interference, we recommend operating the HPLC system separately from the company network.

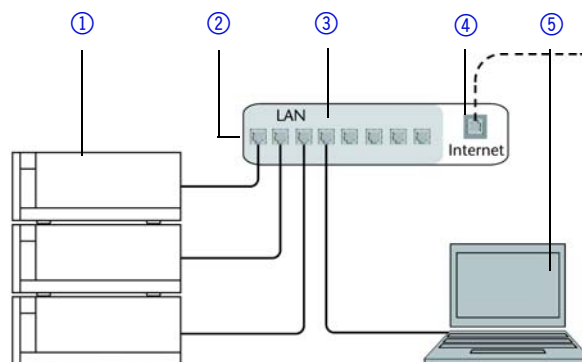You will find patch cables for each device and the router in the accessories kit. To connect the router to a WAN, an additional patch cable is required, which is not supplied within the scope of delivery.

- Prerequisite**
- The computer has been switched off.
 - There is a patch cable for each device and the computer.

- Procedure**
1. Use the patch cable to connect the router and the computer. Repeat this step to connect all devices.
 2. Use the power supply to connect the router to the mains power system.

Configuring the Router

The router is preset at the factory. You will find a label at the bottom side of the router, on which IP address, user name, and

password are printed. These information help to open the router configuration.

- Procedure**
1. To open the router configuration, start your Internet browser and enter the IP address (not for all routers).
 2. Enter user name and password.
 3. Configure the router as DHCP server.
 4. In the router configuration, check the IP address range and make changes if necessary.

Result Once the router has assigned IP addresses to all devices, the chromatography software can be used to remotely control the system.

Integrating the LAN into a Company Network

A network administrator can integrate the LAN into your company network. In this case you use the WAN port of the router.

Prerequisite There is a patch cable for the connection.

- Procedure**
1. Check that the IP address range of the router and of the company network do not overlap.
 2. In case of an overlap, change the IP address range of the router.
 3. Use the patch cable to connect the router WAN port to the company network.
 4. Restart all devices, including the computer.

Controlling Several Systems Separately in a LAN

Devices connected to a LAN communicate through ports, which are part of the IP address. If more than one HPLC system is connected to the same LAN and you plan on controlling them separately, you can use different ports to avoid interference. Therefore, the port number for each device must be changed and this same number must be entered into the device configuration of the chromatography software. We recommend to use the same port number for all devices in the same system.

Note: The port is set to 10001 at the factory. You must use the same numbers in the device configuration of the chromatography software as in the device, otherwise the connection fails.

- Procedure**
1. Find out port number and change it on the device.
 2. Enter the port number in the chromatography software.

Result The connection is established.

Changing the Swing Direction of the Door

The thermostat is delivered with a door which is opened on the right-hand side. Optionally, the swing direction of the door can be changed. To do so, please contact the Technical Support.

Installing the Cartridge for Precolumn Tempering


For precolumn tempering, a cartridge can be ordered separately. With the help of the cartridge, the solvent temperature can be adjusted to the thermostat temperature before it enters the column. This way, a temperature gradient on the column can be avoided and the reproducibility and separating performance are increased.

The cartridge for the precolumn tempering is fixed to the ceiling of the thermostat inside of the thermostat with two screws.

- Prerequisite**
- The device has been switched off.
 - The power plug has been pulled.

Tool Allen wrench, size 2.5

Procedure

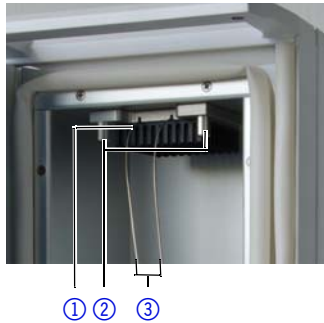Process	Figure
<ol style="list-style-type: none"> 1. Tighten the cartridge ① with two screws ②. 2. Connect the autosampler to one of the two capillaries ③ of the precolumn tempering. 3. Connect the column to the other capillary. 	 <p>Fig. 4 Precolumn cartridge</p>

Installing the Cartridge of the Precolumn Tempering

WARNING! Danger of burns! 185 °F (85 °C) hot surfaces in the column compartment. Wear protective gloves during work in the column compartment.

The separating column is connected directly to the precolumn tempering.

Note: Only one separating column can be connected to the precolumn tempering.

Legend

- ① Precolumn tempering
- ② Separating column

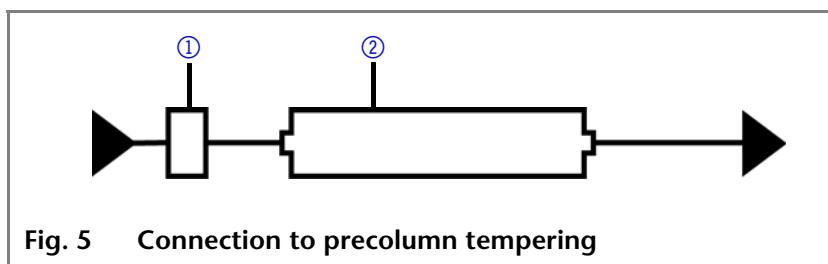
Fig. 5 Connection to precolumn tempering

Operation

In this chapter, you find information which is relevant for operating the thermostat.

Initial Start-Up

Use this checklist to determine whether the detector is ready for the initial start-up:

- Device is positioned in the correct location.
- The power plug has been pulled.

If the device is part of an HPLC system, you should also note the following:

- The network connection to the router is established.
- The chromatography software has been installed by KNAUER or a company authorized by KNAUER.

Switch-On

Prerequisite The installation has been completed.

WARNING! Beware of explosion! High concentration of organic solvent vapor may lead to explosion.

Can cause serious injury. Do not allow flammable chemicals to leak.

CAUTION!

Possible damage to the device caused by condensed water!

Allow device to acclimate for 3 h, before it is connected to power supply and taken into operation.

- Procedure**
1. Connect the power supply of the device to the mains.
 2. Switch the power supply on.
 3. Enter the temperature via chromatography software or the operating panel.

Stand-by If the thermostat is not operating, the heating and cooling elements are switched off. After around 4-5 minutes after operation has ended, the fan inside of the thermostat is switched off as well. Only the outer fan and the controller continue to consume energy.

Controlling with Chromatography Software

The device has to be connected to a computer via LAN to be able to control it with the software packages OpenLAB EZChrom Edition or ClarityChrom.

You find a detailed description on the chromatography software in its accompanying user manual.

Controlling with Operating Panel

Alternatively, the device can be controlled with an optionally available remote operating panel. To control the device using the remote operating panel, connect the operating panel with the keyboard connector on the rear of the device.

Legend

- ① Display
- ② Buttons

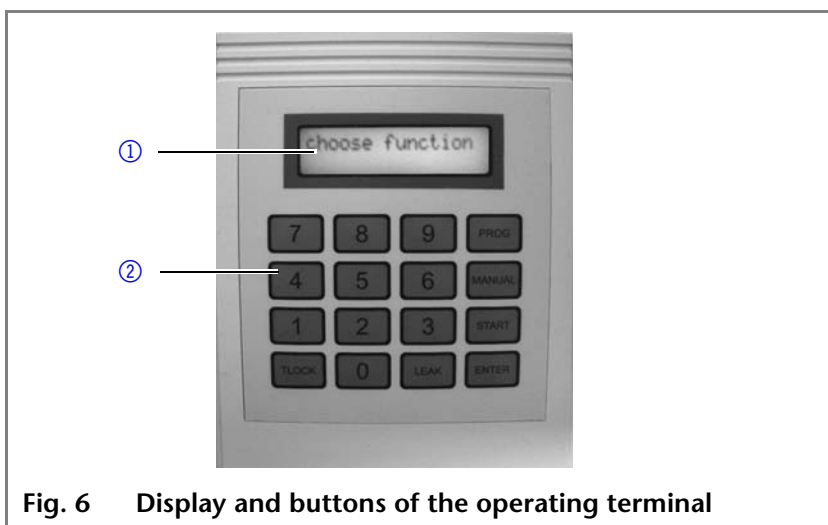
Fig. 6 Display and buttons of the operating terminal

The operating panel consists of a 4x4 matrix keyboard and a double-spaced backlit display.

Intensity of the display You can adjust the intensity of the display through a small opening on the bottom of the operating panel. Use a screwdriver to adjust the intensity.

Stand-by If "choose function" is displayed, the device is in stand-by.

Self-Testing

After starting-up the thermostat, the device name, firmware version and serial number are displayed on the operating panel. The device performs a self-test which checks the individual components of the device. If the test has been completed successfully, "choose function" is displayed.


Remote Operation Using Chromatography Software

'Remote' operation If the chromatography devices are operated with chromatography software, they cannot be controlled with the operating panel. During remote operation using chromatography software, *Remote* is displayed in the status bar.

Leak Sensor Settings

Upon delivery, the default sensitivity of the leak sensor is set to 10 % (operating panel display) or *low* (software display).

Procedure

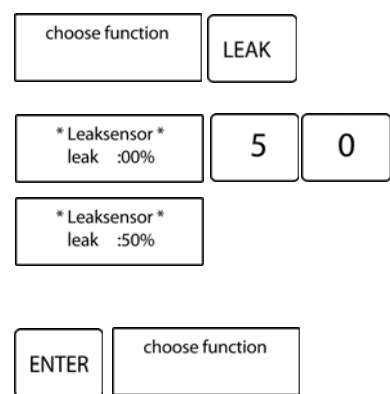
Process	Figure
<ol style="list-style-type: none"> 1. When <i>choose function</i> appears, press the <i>LEAK</i> key. 2. In the <i>Leaksensor</i> menu, use the numeric keys to enter a solvent sensitivity within the range of 00-99 %, e.g. 50 %. 3. Press <i>ENTER</i> to save settings. 	 <p>Fig. 7 Leak sensor settings</p>

Note: If the leak sensor sensitivity is set to 00 %, the leak sensor is switched off.

Setting the T-Lock Temperature

In this menu, you can define an upper temperature limit between 5 °C and 85 °C. If this limit is reached, the device interrupts the power supply of the heating and cooling elements and displays a message. This is to protect your column from accidental overheating.

Procedure

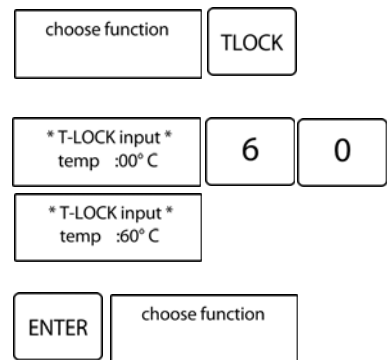
Process	Figure
<ol style="list-style-type: none"> When <i>choose function</i> appears, press the <i>TLOCK</i> key. In the <i>T-Lock input</i> menu, use the numeric keys to enter a temperature within the range of 05-99 °C, e.g. 60 °C. Press <i>ENTER</i> to save settings. 	 <p>Fig. 8 Changing the T-Lock temperature</p>

Note: Entering values < 05 °C or > 85 °C deactivates the T-Lock functionality.

Entering a Fixed Temperature

If the thermostat shall be operated at a fixed temperature, it can be set as follows:

Procedure

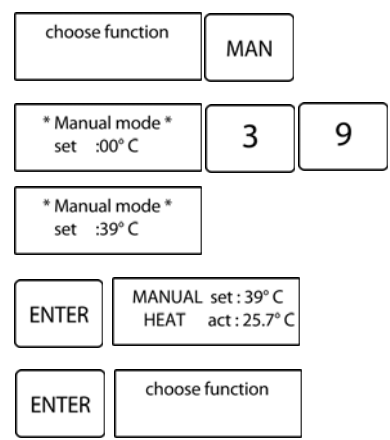Process	Figure
<ol style="list-style-type: none"> When <i>choose function</i> appears, press the <i>MAN</i> key. In <i>Manual mode</i>, press the numeric keys to enter the temperature, e.g. 39 °C. Press <i>ENTER</i> to heat the column thermostat. Press <i>ENTER</i> to end manual temperature settings and to set the device to standby. 	 <p>Fig. 9 Setting the temperature</p>

Note: The current temperature is displayed under "act". The status of the device is designated as follows:

- HEAT: The device is in heating mode.
- COOL: The device is in cooling mode.
- IDLE: The device does not heat and not cool.

Changing the Temperature

If a fixed temperature has been set, it can be adjusted later on during operation.

Procedure

Process	Figure
<div>1. Reducing the temperature: When <i>MANUAL</i> appears, press the <i>TLOCK</i> key to reduce the temperature by 1 °C. Keep <i>TLOCK</i> pressed to further reduce the temperature incrementally.</div> <div>2. Increasing the temperature: When <i>MANUAL HEAT</i> appears, press the <i>LEAK</i> key to increase the temperature by 1 °C. Keep <i>LEAK</i> pressed to further increase the temperature incrementally.</div>	<div><div><div>MANUAL set : 39° C HEAT act : 25.7° C</div><div>TLOCK</div></div><div><div>MANUAL set : 38° C HEAT act : 25.7° C</div><div>TLOCK</div></div></div> <div><div><div>MANUAL set : 39° C HEAT act : 25.7° C</div><div>LEAK</div></div><div><div>MANUAL set : 40° C HEAT act : 25.7° C</div><div>LEAK</div></div></div> <div><div>Fig. 10</div><div>Changing the temperature</div></div>

Entering Timed Temperature Programs

With the programming functionality, it is possible to set time-limited temperature segments and up to 99 programs.

Procedure

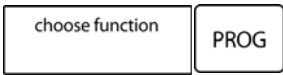
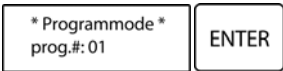
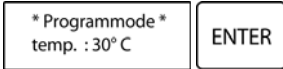
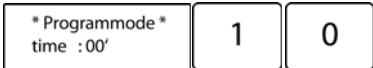
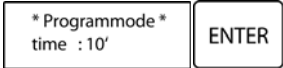
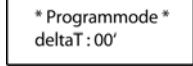
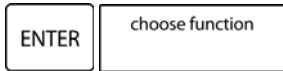
Process	Figure
1. When <i>choose function</i> appears, press the <i>PROG</i> key.	
2. In the <i>Programmode: progmenu</i> , use the numeric keys to enter a number for the program, e.g. 01.	
3. Confirm the settings with <i>ENTER</i> .	
4. In the <i>Programmode: temp menu</i> , enter the temperature in °C (05-85 °C, e.g. 30) eingeben.	
5. Confirm the settings with <i>ENTER</i> .	
6. In the <i>Programmode: time menu</i> , enter the program duration (01-98 min., e.g. 10).	
7. Confirm the settings with <i>ENTER</i> .	
8. In the <i>Programmode: deltaT menu</i> , no changes have to be made.	
9. Finish settings with <i>ENTER</i> . Simultaneously, this saves the program.	

Fig. 11 **Creating a temperature program**

Note: To delete old programs, it is sufficient to overwrite them.

Linked Programs

If multiple programs are set, the device automatically runs through all programs subsequently.

Note: Remember to allow enough time for the separate temperature increments when setting the programs. This ensures that the target temperature can be reached. Allow for a heating rate of 2-3 °C/min. The heating rate is reached when the recommended ambient temperature (17–28 °C) is met. If the selected time segment for a temperature increment is too short, the set temperature for said increment may not be reached. This is due

	to the system continuing with the subsequent temperature increment after the set time has elapsed.
Changing the program number	If you want to change the program number of multiple runs manually during a run, use the <i>TLOCK</i> key to move to the prior key and the <i>LEAK</i> key for the subsequent program number.
Loop functionality	To combine temperature gradients to a loop, finish with the program time = 00. Program # 00 is switched back to and the sequence is repeated. To use the Loop functionality, the first program of the sequence must be program # 00. Note that time = 00 is a default setting of every program. If the settings have not been changed after the last program, the loop functionality is activated automatically.
Automatic sequence end	If not sequence end has been set, the CT 2.1 runs through all entered programs and then jumps to program #00 in loop mode. If only one part of all set programs should be run through or no loop be performed, add an additional program with temp = 00 and time = 1. By this, the sequence is stopped and the following programs are not run through.
Infinite functionality	If you want to continue the last program infinitely, enter time = 99 in the last program to be run through. By this, the temperature of the last run is maintained indefinitely.
Manual sequence end	With the <i>ENTER</i> key, you can exit the program functionality at any time.

Programming a Temperature Gradient

The column thermostat CT 2.1 can connect single temperature plateaus and calculate defined temperature gradients. For this, a heating and cooling period prior to a temperature plateau is calculated using the *deltaT* functionality. By this, the temperature between two plateaus is aimed at, not with a maximum heating or cooling rate, but defined via the rise time or the decay time.

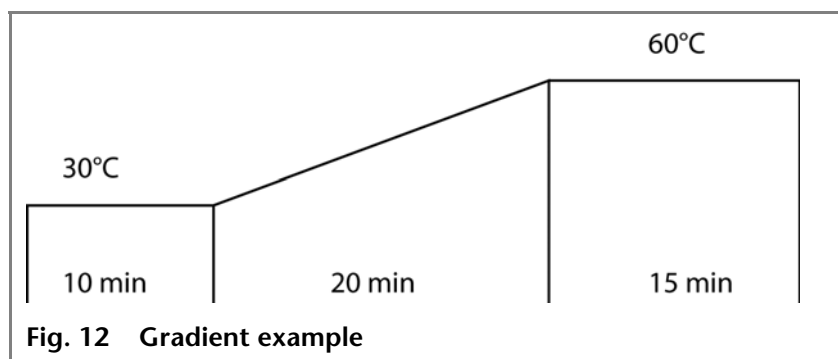
Fig. 12 Gradient example

To reproduce the programming of the gradient pictured above, proceed as follows:

Procedure

Process	Figure
<p>Set the program # 00 with 10 min at 30 °C as described in "Entering Timed Temperature Programs". Afterwards, perform the following settings for program # 01:</p> <ol style="list-style-type: none">1. When <i>choose function</i> appears, press the <i>PROG</i> key.2. In the <i>Programmode: prog</i> menu, use the numeric keys to enter a number for the program, e.g. 01.3. Confirm the settings with <i>ENTER</i>.4. In the <i>Programmode temp</i> menu, enter the temperature 60 °C.5. Confirm the settings with <i>ENTER</i>.6. In the <i>Programmode time</i> menu, enter the program time of 15 min.7. Confirm the settings with <i>ENTER</i>.8. In the <i>Programmode deltaT</i> menu, enter 20 min. as the time span allocated to reaching the temperature plateau.9. Finish settings with <i>ENTER</i>. Simultaneously, this saves the program.	<p>Fig. 13 Programming a gradient</p>


Note: To exit the gradient functionality during a run and skip to the next temperature plateau, press *ENTER*.

Starting a Program or a Sequence

To start a program or a sequence, proceed as follows:

Procedure	Process	Figure
	<div>1. When <i>choose function</i> appears, press the <i>START</i> key.</div> <div>2. In the <i>select programm</i> menu, enter the program by pressing the numeric keys, e.g. 01.</div> <div>3. Press <i>ENTER</i>. The program starts.</div>	<div><div>choose function</div><div>START</div></div> <div><div>select program prog.#: 00</div><div>0</div><div>1</div></div> <div><div>select program prog.#: 01</div><div>ENTER</div></div> <div>Fig. 14 Starting a program</div>

The following parameters are displayed during a run.

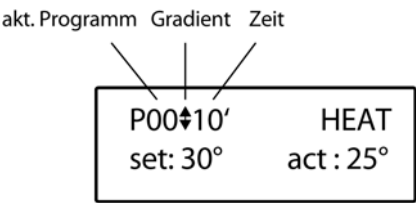P00 current program

Time remaining time in this section of the program

↑ increasing temperature gradient

↓ decreasing temperature gradient

⬮ temperature plateau

settarget temperature

actcurrent temperature

HEAT the thermostat is heating

COOL the thermostat is cooling

Note: The temperature that is displayed may deviate from the target temperature because this closed-loop control is very sensitive.

Functionality Tests

Installation Qualification (IQ)	<p>The customer may request the Installation Qualification, which is free of charge. In case of a request, the Technical Support of KNAUER or from a provider authorized by KNAUER performs this functionality test during the installation. The Installation Qualification is a standardized document that comes as part of the delivery and includes the following:</p> <ul style="list-style-type: none"> ▪ confirmation of flawless condition at delivery ▪ check if the delivery is complete ▪ certification on the functionality of the device
Operation Qualification (OQ)	<p>The Operation Qualification includes an extensive functionality test and must be purchased from the manufacturer. Contact the KNAUER Sales Department to request an offer. The Operation Qualification is a standardized KNAUER document and includes the following:</p> <ul style="list-style-type: none"> ▪ definition of customer requirements and acceptance terms ▪ documentation on device specifications ▪ device functionality check at installation site
Test Intervals	<p>To make sure that the device operates within the specified range, you should test the device using the Operation Qualification at following intervals:</p> <ul style="list-style-type: none"> ▪ Every 3 months: average useful life of more than 5 days/week or 24 hours/day; when operating with buffer solutions or other salt solutions: ▪ Every 6 months: average useful life of 1 to 5 days/week
Execution	<p>The test can be carried out either by the Technical Support of KNAUER or from a provider authorized by KNAUER.</p>

Troubleshooting

- First measures**
1. Check all cabling.
 2. Check device for leaks.
 3. Pay attention to system messages.

LAN

Go through the following steps, in case no connection between the computer and the devices can be established. Check after each step if the problem is solved. If the problem cannot be located, call the Technical Support.

1. Check the status of the LAN connection in the Windows task bar:
 -  Connected
 -  Connection not established

If no connection was established, test the following:

- Is the router switched on?
 - Is the patch cable connected correctly to the router and the computer?
2. Check the router settings:
 - Is the router set to DHCP server?
 - Is the IP address range sufficient for all the connected devices?
 3. Check all connections:
 - Are the patch cable connected to the LAN ports and not the WAN port?
 - Are all cable connections between devices and router correct?
 - Are the cables plugged in tightly?
 4. If the router is integrated into a company network, pull out the patch cable from the WAN port.
 - Can the devices communicate with the computer, even though the router is disconnected from the company network?
 5. In case you own a Control Unit, check the settings in the menu *Setup > Network*.
 - Is *LAN-DHCP* set for controlling?
 - Did the device receive an IP address?
 6. Turn off all devices, router, and computer. Firstly, turn on the router and secondly turn on the devices and the computer.
 - Has this been successful?
 7. Replace the patch cable to the device with that no connection could be established.
 - Has this been successful?
 8. Make sure that the IP port of the device matches the port in the chromatography software.

Possible Problems and Rectifications

Error	Solution
Device cannot be switched on	Inspect the power cable to ensure that it is plugged into the power supply.
The alarm does not sound when a leak occurs	Leak sensor settings (see chapter "Leak Sensor Settings" and "Adjusting the Leak Sensor")
The alarm sounds although no leak can be detected	The leak sensor sensitivity is set too high Leak sensor settings (see chapter "Leak Sensor Settings" and "Adjusting the Leak Sensor")

The target temperature has not been reached	<ul style="list-style-type: none"> ▪ Check if the door is fully closed ▪ Check if the ventilation slots on the back and on the side are unblocked ▪ Check if the fans inside and outside of the device are working ▪ Set a longer plateau phase ▪ Correcting the Temperature (see chapter "Correcting the Temperature")
System error	Switch off the device to reset the memory. Restart.

Error messages

If an error occurs, a signal sounds repeatedly. If error messages are displayed which are not listed below, restart the device once. If error messages are repeated, contact the Technical Support.

Error Message	Measures
ERROR temp sensor	Device error during self-test. Contact Technical Support.
T-LOCK ALARM	Temperature limit reached. Check the set temperature and change, if necessary. If no fault can be found in the settings, contact the Technical Support.
LEAKSENSOR ALARM	Solvent leak. Check capillary connectors at the column.
T-SENSOR ALARM	Device error. Contact Technical Support.
ERROR EEPROM	Device error. Contact Technical Support.

Next steps After troubleshooting is finished, press the *ENTER* key to continue.

Maintenance and Care

Organic solvents are toxic above a certain concentration. Ensure that work areas are always well-ventilated! When performing maintenance tasks on the device, always wear safety glasses with side protection, protective gloves, and a lab coat.

CAUTION!

Performing maintenance tasks on a switched on device can cause damage to the device. Switch off the device and pull the power plug.

Contact the Technical Support

If you have any technical questions regarding KNAUER hardware or software, please use one of the contact options below:

Technical Support hotline:

European hotline	Languages:	Available by telephone
	in German or English:	8 am to 5 pm (CET)
	Phone:	+49 30 809727-111
	Telefax:	+49 30 8015010
E-mail contact:		support@knauer.net

Maintenance Contract

The following maintenance work on the device may only be performed by KNAUER or a company authorized by KNAUER and is covered by a separate maintenance contract:

- Opening the device or removing housing parts

Cleaning and Care

CAUTION!

Intruding liquids can cause damage to the device.

- Place solvent bottles next to the device or in a solvent tray.
- Moisten the cleaning cloth only slightly.

Clean all smooth surfaces of the device with a mild, commercially available cleaning agent or with isopropanol.

Display Clean the display of the operating panel with isopropanol and dry with a smooth, lint-free cloth.

Storage

The device can be stored under the following ambient conditions:

- Temperature range 10-40 °C (50-104 °F)
- Humidity below 90 %, non-condensing

Disconnecting from Power Supply

Prerequisite The device has been switched off.

Procedure

Process

1. Pull the power plug from the socket first and then from the device.
2. Store the power cable together with the device.

Next steps Disconnect the remaining electrical connections from the device and pack them with the device for transport or storage.

Removing a Leak

Prerequisite If liquid enters the device, switch off the device.

Auxiliary material Cloth for drying the leak sensor

Procedure

Process

1. Remove the leak
2. Dry the leak tray.
3. Confirm the error message via the operating panel or the chromatography software.

Next steps Put the device into operation.

Correcting the Temperature

If the temperature displayed on the CT 2.1 deviates from the actual temperature, you may correct the temperature yourself.<

Temperature - Auto-Calibration

This functionality allows fully automated calibration of the device. The procedure takes 18 hours and ends with a verification. During auto-calibration, the device runs through the following temperature values:

85 °C, 70 °C, 60 °C, 50 °C, 40 °C, 30 °C, 20 °C, 10 °C, 05 °C

Note: The room temperature for auto-calibration has to be stable at 22-23 °C.

Procedure

Process

1. Press *TLOCK*.
2. Enter 0 °C and finish with *ENTER*.
3.

9

5

6

 and finish with

0

 .
4. If the device displays *choose function* again, switch off the device. Restart after 15 seconds to save the changes.

Correcting the Temperature Manually

Note: The verification should be performed by the Technical Support only, in order to ensure the precision of the device. If, however, you want to correct the temperature with a verified measuring instrument, proceed as follows:

The temperature can be corrected manually at the following points:

05 °C, 10 °C, 20 °C, 30 °C, 40 °C, 50 °C, 60 °C, 70 °C

Adjust the temperatures subsequently. Write down the difference between the current temperatures displayed on your measuring instrument and the temperatures displayed on the CT 2.1. The difference is your correction value.

Example If the temperature set in the CT 2.1 is 5 °C and the temperature on your measuring instrument is 5,1 °C, the correction value is -0,1.

Procedure	Process
	<p>1. Enter the combination in the menu <i>choose function</i></p> <div data-bbox="628 931 916 1001" style="border: 1px solid black; padding: 5px; display: inline-block;"> <div style="border: 1px solid black; padding: 2px 10px; margin: 0 5px;">9</div> <div style="border: 1px solid black; padding: 2px 10px; margin: 0 5px;">4</div> <div style="border: 1px solid black; padding: 2px 10px; margin: 0 5px;">9</div> </div> <p>The menu <i>set value corr.</i> is displayed with the first temperature value. If the temperature value does not have to be changed, switch to the next value with <i>LEAK</i>.</p> <p>To change the temperature, press <i>ENTER</i> and enter the correction value. Change the algebraic sign to minus with <i>TLOCK</i> or to plus with <i>LEAK</i>.</p> <p>2. Close every entry <i>ENTER</i>.</p> <p>3. All values can be checked with <i>TLOCK</i> (scrolls down) and <i>LEAK</i> (scrolls up).</p> <p>4. To leave the menu, press <div data-bbox="960 1384 1051 1453" style="border: 1px solid black; padding: 2px 10px; display: inline-block;">1</div>.</p> <p>5. Switch off the thermostat for one minute to activate the changes.</p>

Note: Maximum deviations that can be entered are +/-2,5 °C.

Adjusting the Leak Sensor

The leak sensor functionality should be checked annually. This is possible with the optional operating panel only.

To test the leak sensor functionality, proceed as follows:

Necessary materials Syringe; methanol: 25 µl and 100 µl, cellulose

State of the device The device is in an upright position, its door closed, and running with a temperature of between 20 °C and 30 °C.

Methanol injection Spray the methanol onto a piece of cellulose and immediately place it inside of the thermostat. Close the door.

Test results 25 µl methanol - no alarm

100 µl methanol - alarm is activated

Note: Between both tests, open the door for a few minutes to allow the methanol used in the first test to dissipate.

Correcting the leak sensor If the test does not yield the expected result, enter 957 into the operating panel. Using a potentiometer, set the value "L" on the rear of the device.

Legend

① Potentiometer

Fig. 15 Potentiometer

Technical Data

Thermostat	Tempering	microprocessor-controlled Peltier device for heating and cooling, the fan supports two-way air circulation.
	Temperature range	5–85 °C
	Heating/cooling rate	2 °C/min
	Temperature precision	+/- 0,2 °C
	Temperature stability	+/- 0,1 °C
Column compartment	Pre-heating of the mobile phase	passive tempering with an optional cartridge
	Number of columns	up to 4
	Column size	analytical and preparative columns with up to 20 mm ID and 350 mm length.
	Dimensions of the column compartment	90 x 390 x 47 mm (W x H x D)

Communication	Control	<ul style="list-style-type: none">▪ LAN▪ optional Control Unit	
	Programming	temperature and gradients	
Technical parameters	Housing door	variable swing direction with door sensor	
	GLP	via software: serial number, firmware version, date of installtion, date of last service, main-board run time	
	Safety	self-test and auto-calibration upon switch-on; adjustable switch-off temperature	
	Operating panel	optional	
	Leak sensor	gas sensor with adjustable sensitivity and acoustic signal	
	Degree of protection	IP-20	
	Dimensions	150 mm x 470 mm x 310 mm (Width × Height × Depth)	
	Weight	8.4 kg	
	Power supply	Power input	90–230 V
		Output	50–60 Hz
Maximum power consumption		100 Watt	
Operating Environment	Recommended ambient temperature	17–28 °C (62.6–82.4 °F)	
	Humidity	below 90 %, non-condensing	
	Height above sea level	maximum 2000 meters	

Repeat Orders

Find further information regarding spare parts and accessories at www.knauer.net.

Modules

Name	Order number
Column Thermostat CT 2.1	A05852

Accessories and Spare Parts

Name	Order number
Cartridge ID 0,25 mm	A05852-2
Cartridge ID 0,1 mm	A05852-3
Operating panel	A05852-1
Accessory kit	F05852
Installation Qualification EN	VIQC21
Operation Qualification EN	VOQC21

Qualification
documents

Disposal

Hand in old devices or disassembled old components at a certified waste facility, where they will be disposed of properly.

AVV Marking in Germany

According to the German "Abfallverzeichnisverordnung" (AVV) (January, 2001), old devices manufactured by KNAUER are marked as waste electrical and electronic equipment: 160214.

WEEE Registration

KNAUER as a company is registered by the WEEE number DE 34642789 in the German "Elektroaltgeräteregister" (EAR). The number belongs to category 8 and 9, which, among others, comprise laboratory equipment.

All distributors and importers are responsible for the disposal of old devices, as defined by the WEEE directive. End-users can send their old devices manufactured by KNAUER back to the distributor, the importer, or the company free of charge, but would be charged for the disposal.

Solvents and Other Operating Materials

All solvents and other operating materials must be collected separately and disposed of properly.

All wetted components of a device, e. g. flow cells of detectors or pump heads and pressure sensors for pumps, have to be flushed first with isopropanol and then with water before being maintained, disassembled or disposed.

Legal Information

Warranty Conditions

The factory warranty for the device is valid for 12 months after the date of dispatch. All warranty claims shall expire in the event that any unauthorized changes are made to the device.

During the warranty period, any components with material or design-related defects will be replaced or repaired by the manufacturer free of charge.

This warranty excludes the following:

- accidental or willful damage
- damage or errors caused by third parties that are not contractually related to the manufacturer at the time the damage occurs
- wear parts, fuses, glass parts, columns, light sources, cuvettes and other optical components
- damage caused by negligence or improper operation of the device and damage caused by clogged capillary
- packaging and transport damage

In the event of device malfunctions, directly contact the manufacturer.

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Transport Damage

The packaging of our devices provides the best possible protection against transport damage. Check the devices for signs of transport damage. In case you notice any damage, contact the Technical Support and the forwarder company within three workdays.

HPLC Glossary

The following table gives you an overview of abbreviations and terminology used in this manual.

Term	Definition
capillary	Thin metal or PEEK tube that connects components or devices in an HPLC system.
column	The column is filled with a matrix. This matrix can separate individual substances by interacting differently with the analytes in a solvent.
chromatography	Procedure for the separation of a mixture of substances with the help of a mobile and a stationary phase.
GLP	A quality control mechanism for laboratories (Good Laboratory Practice).
HPLC	High Pressure Liquid Chromatography
IP address	Unique address of a sender or receiver in local networks or on the internet (Internet protocol).
IQ	Proof of adequate installation of a device.
LAN	Local area network (LAN). Network in which HPLC devices and a computer are connected with each other to access the devices.
LAN-DHCP	Local computer network with automatic configuration (Local Area Network - Dynamic Host Configuration Protocol), consisting of an ethernet card and a LAN connection cable (patch cable):
mobile phase	The flowing agent that transports the substances to be separated or isolated through the column.
OQ	Extensive functionality test of the single devices of an HPLC system (Operation Qualification).
Peltier device	Electrochemical converter that generates a temperature difference when current from one side of the converter to the other.
port	Internal computer address of the interface of the local network.
PQ	The certificate of performance of an HPLC system (Performance Qualification) based on a standardized testing environment.

Term	Definition
pump	Device which conveys the solvent through a chromatographic system with a controlled stream.
operating panel	Optional control device for the CT 2.1, suitable for stand-alone use.
pre-column tempering	With this component, the temperature of the mobile phase is adjusted to the temperature of the thermostat before it enters the column.
remote	Device control via chromatography software. The device cannot be controlled locally (via the operating panel).
router	Device in a computer network which checks and forwards data packages.
RSD	Standard deviations in retention time and peak area found during system tests of the HPLC system (Relative Standard Deviation).
solvent	see mobile phase.
stationary phase	The matrix in the separating column of a chromatographic system against which the substance mixture separates.

Index

A

Accessories 7
Ambient conditions 11, 29
AVV marking 34

C

Capillary 36
Care 29
Chromatography 36
Chromatography software
 control 18
 remote control 19
ClarityChrom® 18
Cleaning 29
Column 36
Column thermostat
 operating panel 18
 set temperature 20, 25
Contact data 29
Control
 chromatography software 18
 operating panel 18
 remote 37

D

declaration of conformity 40
decontamination 9
Degree of protection 33
Dimensions 33
Disposal 34
Door
 swing direction 16

G

GLP 33, 36

H

HPLC 36

I

Initial start-up 17
Intended Use 5
IQ 26

L

LAN 12, 36
 port 15
 problems 26
 router 14
 settings 13
 setup 14
 troubleshooting 26
leak 8
Leak sensor 30
Leak sensor settings 19

M

Maintenance
 maintenance contract 29

O

OpenLAB EZChrom 18
Operating panel 37
 column thermostat 18
Operation
 Remote 19
OQ 26, 36

P

Peltier device 36
port (LAN) 15
power
 strip 8
 supply 8
Power consumption 33
power supply
 cable 8
PQ 36
Precolumn tempering 16, 37
Problems and rectifications 27
professional group 7
Programming a temperature gradient 23
Pump 37

R

Remote operation 19
 chromatography software 19
Router 37
 router (LAN) 14
RSD 37

S

safety equipment 7

- Safety for Users **9**
- Self-test **18**
- Service **29**
- Set the temperature **20**
- Setting the temperature **25**
- Solvent **37**
- solvent
 - flammability **8**
 - line **8**
 - self-ignition point **8**
 - tray **8**
- Spare parts **7**
- Stationary phase **37**
- Storage **29**
- Swing direction of the door **16**
- Switch-on **17**

T

- Technical Data **32**
- Technical data
 - weight **33**
- Technical Support **29**
- test
 - Installation Qualification **26**
 - Operation Qualification **26**
- transport damage **35**
- Troubleshooting **26**
 - LAN **26**

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Model/Type Reference Column Thermostat CT 2.1 Product no. A05852

The product complies with the following standards:

Machinery	Machinery Directive 2006/42/EC IEC 60799:1999
EMC	EMC Directive 2004/108/EC IEC 61000-3-2:2006 + A1 + A2 IEC 61326-1:2013
Disposal	RoHS Directive 2002/95/EU WEEE Directive 2002/96/EU
Safety	Low Voltage Directive 2006/95/EC IEC 61010-1:2011 IEC 61010-2-010:2003 IEC 61010-2-081:2001 + A1

The product was tested with a typical configuration. The mark of conformity has been applied to the rear panel.

Date Berlin, 01.10.2014

Alexandra Knauer (CEO and owner)

© KNAUER Wissenschaftliche Geräte GmbH
All rights reserved.
The information in this document is subject to change
without prior notice. Translation of the original German
edition of this manual, version 1.1
Last manual update: 2015-07-28
Printed in Germany on environmentally friendly paper
from sustainable forests.

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® AZURA is a registered trademark of